

# Crystal Eyes



## Plug-and-play crystallisation monitor

### Features and Benefits

- Monitor and log turbidity, temperature and pH
- Probes for any size vessel
- Single compact portable unit



A complete  
portable system for  
crystallisation studies

For more information or a demonstration please contact us  
[www.helgroup.com](http://www.helgroup.com) T: +44 (0) 20 8441 6778 E: [marketing@helgroup.com](mailto:marketing@helgroup.com)

Easily determines the solubility and super-solubility temperatures as you heat and cool your solution

### Crystal EYES includes

- Turbidity and temperature probes chosen specifically for your application
- Electronics and PC with data handling software
- Optional pH monitoring

### Ease of use

- Insert the probes into your vessel
- Heat and cool your solution
- Data is automatically displayed and recorded

### Choose

- The probe required size for your application
- The material best equipped for your chemistry and budget

### Data handling

- Real time data display and capture
- The data is automatically saved to a data file
- The data is readily exported to Excel for analysis




Probes for vessels as small as test tubes and as large as bench scale reactors

#### Sizes:

Test tubes (10cm) Small Vessels (15cm)  
Bench Scale (50cm)

#### Materials:

Peek, Stainless steel Hastelloy

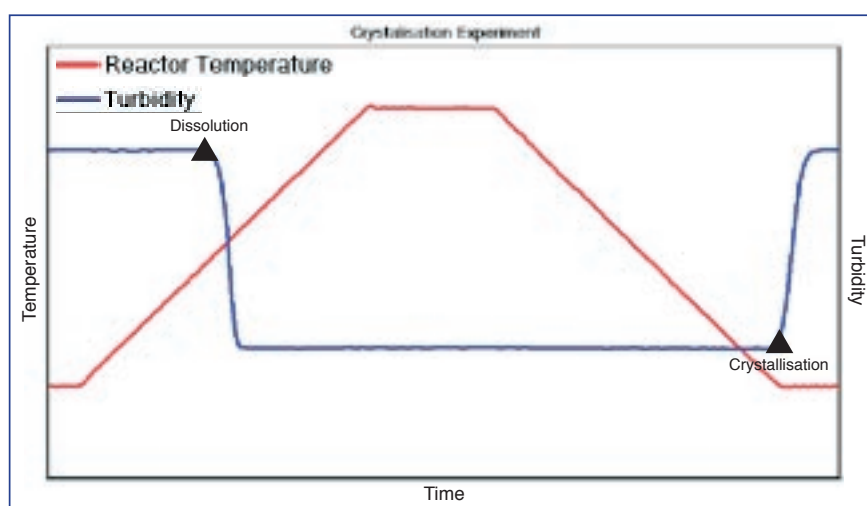

Typical trace data from a Crystal Eyes solubility experiment

For more information or a demonstration please contact us  
[www.helgroup.com](http://www.helgroup.com) T: +44 (0) 20 8441 6778 E: [marketing@helgroup.com](mailto:marketing@helgroup.com)